

Bite Stick Technique

Esthetic rehabilitation requires detailed planning and good communication. Making sure that your patient's facial midline and the dental midline coincide and determining the cant and horizontal incisal edge position of the central incisors, has always been a challenge to communicate, until now!


The technique requires the use of a flexible bite material such as Blue Mousse™ and a swabbing stick. Or “Bite Stick. Or (**Symmetry Facial Plane Relator, Clinicians Choice New Milford, CT. 800-265-3444**)

With the teeth closed the bite material holds the stick in position on the lower incisors as it is aligned to the desired horizontal incisal plane. Use the patients facial features to create the final result. If you need to indicate the mid-line it can be scribed on the stick. Make sure enough material is used to hold the stick in place on the lower incisal edge. This step is not intended to be used as a bite for occlusion purposes. A separate bite should be included for occlusion

Use this easy technique and increase productivity and predictability.

Advantages/Benefits

- No more slanted incisal planes
- Provides technical information for predictable results.
- Your technician will love it.

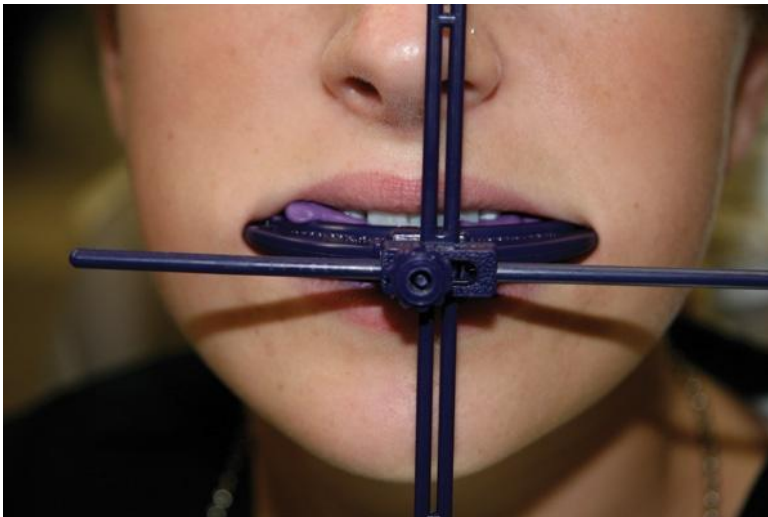

PRO-Craft
DENTAL LABORATORY, INC
INSTITUTE OF DENTAL TECHNOLOGY

25791 Jefferson Ave. Murrieta, CA. 92562 www.pro-craft.com 877-484-3522